


CAMERON LIFT CHAIR

The Cameron lift chair makes life easy with a powered lift function for those with mobility difficulties, plus powered recline, a heated seat backrest and handy USB outlet





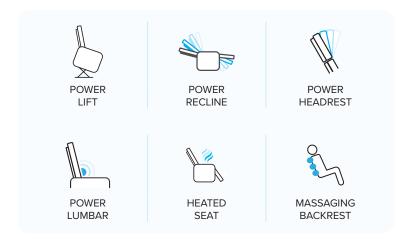





CAMERON LIFT CHAIR

UPHOLSTERY OPTIONS

STOCKED IN:



ZEUS FABRIC (BALTIC)



scan QR code to view this product on synargy.com

FEATURES & BENEFITS		
FRAME	FSC certified Pine and Plywood frame	A sturdy frame constructed with sustainably sourced timber provides a solid and long- lasting base to the suite
FOUNDATION	Heavy gauge steel Zig Zag springs	The foundation delivers long-lasting support & the first layer of comfort
FOAM	Premium high-density D30 foam topped with a polypropylene fibre layer	Layered high & low density foam delivers seating that lasts without compromising on comfort
MECHANISM	Tested up to 150kgs at 25,000 cycles	Warrantied for use up to 150kgs, this robust recliner mechanism has been crafted to last that's 30kg more weight-loading than standard mechanisms!
POWERED LIFT	Powered lift function	Powered lift allows for easy exit out of the chair - perfect for those with mobility difficulties
TRIPLE POWERED	Powered recliner, headrest and lumbar	Recliners, headrests and lumbar are adjustable at the touch of a button - the ultimate for customised comfort!
AIR MASSAGE	Airbag massage in seat back	Kick back and relax with a soothing 30-minute air-powered massage in the seat back. Allow 10-minutes between massage cycles.
HEATED SEAT	Heated seat backrest	Beat the winter chill with a heated backrest that provides soothing warmth and gentle muscle relief at the press of a button.
USB OUTLET	USB outlet in the remote control	Convenient charging at your fingertips ensures you never need to be separated from your device. Charge smartphones, tablets, and other small gadgets while comfortably seated. USB port located in the remote control
UPHOLSTERY	Available in Zeus fabric	Stocked in a Baltic coloured Zeus fabric











CAMERON LIFT CHAIR

REMOTE CONTROL USER GUIDE Head Rest buttons headrest Used to adjust the angle of the head rest Lift buttons Used to control the power lift function of the chair lift **Back buttons** back Used to adjust the amount of back rest recline Home button lumbar 😭 **Lumbar buttons** Push and hold to restore the chair back to the Used to adjust the position of the lumbar area default or "home" position. Holding the home button past the default position will activate the lift function to assist the user out of the chair. → footrest mode **Footrest buttons** Mode button Used to control the extension of the footrest Used to cycle between different massage modes: · Green Light: Pulse Mode Z O Z O Z O Z O • Blue Light: Wave Mode · Red Light: Alternating n e n n е e 1, 2, 3, 4 buttons 2 3 4 Used to activate or deactive particular massage areas in the backrest: Zone 1: Lower back • Zone 2: Middle lower back • **Power button** Zone 3: Middle upper back • Used to turn the back rest massage on or off. The massage Zone 4: Upper back • will run for 30-minutes before automatically turning off. Allow 10-minutes between massage cycles. Heat button Used to turn back rest heating on or off. The heating will run for 30-minutes before automatically turning off. Allow 10-minutes between heat cycles.







